


global[®] Collect[®]

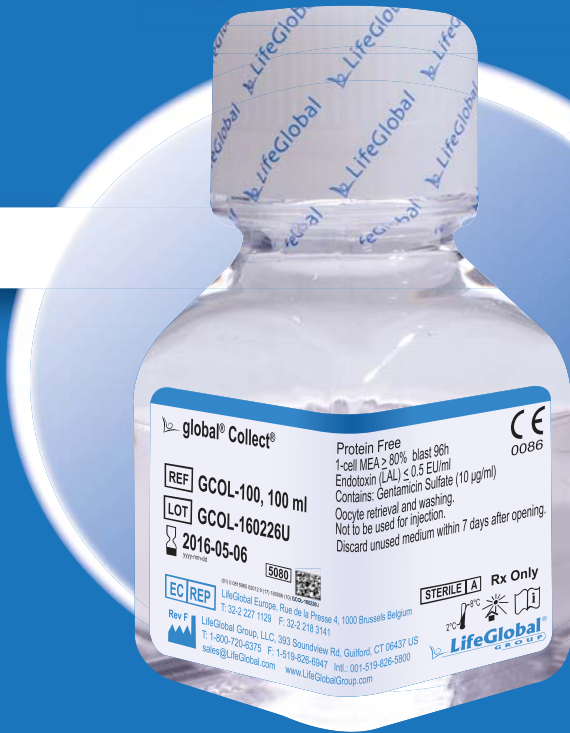

A unique formulation based on global[®]

global[®] Collect[®] is a HEPES-buffered medium based on global[®] formulation for retrieval, holding and washing of cumulus-oocyte complexes after collection to be used at 37°C.

- Contains increased glucose sufficient to satisfy the metabolic requirements of the cumulus cells
- Available in a 500 mL format for high volume users
- Maintain same chemical environment throughout all stages of oocyte handling and culture with global[®] family
- High Quality of Manufacturing
- Stringent Quality Control

global[®] Collect[®] contains the antibiotic gentamicin sulfate. Appropriate precautions should be taken to ensure that the patient is not sensitized to this antibiotic. The pH of global[®] Collect[®] is balanced for use at 37°C. global[®] Collect[®] contains only a low concentration of sodium bicarbonate and should not be gassed, or used for culture under CO₂.

T: 1-800-720-6375 F: 1-519-826-6947 Intl.: 001-519-826-5800
sales@LifeGlobal.com www.LifeGlobalGroup.com

©2017 LifeGlobal® Group, LLC

 **LifeGlobal[®]**
GROUP

global_Collect-2017-08-16